

4IN1 MACHINE

Diode laser E-light ND YAG RF



Real shot display

01



Diode laser handle

02



Elight handle

03



ND YAG Laser handle

04



RF handle

elight ipl 808 nm rf

	Elight808nm RF
	1) 2) 3) 4) 5) 6) 7) 8) 9)

DIODE LASER

755/808/1064nm diode laser, for all type skin hair removal, painless

755 NM

The Power of Diversity



755nm

808 NM

Up To Speed



808nm

1064 NM

Dedicated to Darker Skin



1064nm



One minute cooling

Suitable for hair removal on any part of the body



Arms



Legs



Face



Armpits



Back



Chest




How does it work

E-light IPL for Hair removal/Skin rejuvenation/Vascular Job/Pigmentation treatment/Acne clearance

Face and body with big and small treatment head.

Through the scientific design of diode laser handle, the high-energy laser beam instantly penetrates the skin and reaches the hair follicle during treatment, output light is absorbed by melanin in the hair follicle, so destroying the hair follicle and reach permanent hair removal. The laser pulse width is set shorter than the time that human body's feeling time, a single shot of light is over, but the human body fell painless, then could achieve the painless treatment feeling, to improve the customer's comfortable experience

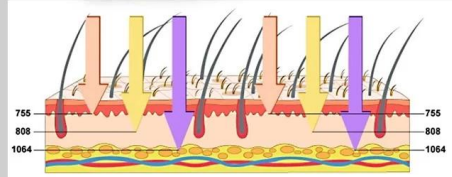

Refuse To Repeat And Leave No Scum



1 week-2 weeks (2 times a week)	3 week-4 weeks (2 times a week)	5 week-8 weeks (2 times a week)
75% hair loss Start to grow slowly Smaller pores	Hair becomes thinner and less Rejuvenate tender skin's smooth body	Select according to hair growth status. Use it selectively, say goodbye Repeated hair removal.

IPL Handle

- ◎ **High-quality IPL**

Handpiece with super cooling system

- ◎ **Optional Spot Size**

15*50mm, 15*35mm, 12*30mm

- ◎ **Long shots time**

IPL: 150,000 times

- ◎ **7 filters**

420/510/560/585/640/690/755nm
Standard: 420/510/560/690nm



420nm	510nm	560nm	585nm	640nm	690nm	755nm
Acne removal Vascular treatment	Pigmentation therapy, Remove freckles	Skin rejuvenation wrinkle removal	Remove Red blood silk	Hair removal (light skin)	Hair removal (brown skin)	Hair removal (black skin)
						



ND YAG

5 Treatment Heads

532: Remove colored tattoos.
 755 honeycomb head: used to remove age spots, freckles, and chloasma. little damage to skin
 Adjustable parameter 1064: for tattoo removal, eyebrow clearing
 Non-adjustable parameter 1064: used for removing tattoos and clearing eyebrows
 1320: large area treatment, large spot, fast treatment speed and good effect



1064nm
(Size Optional)

532nm
(Original)

1320nm
(Blackface doll)

755nm
(Honeycomb)

1064nm
(Original)



ND YAG ADUANTAGE

Imported crystal bonding laser rod, best quality , shot more than 1.5 million times
 Counter on the handle, more easy to get the total shot time
 LED red indicator, more easy to aim on the treatment area.
 5 treatment tips for multi fuctions: tattoo removal, skin rejuvanation, pigment removal, carbon peeling....



RF Handle

Suitable to 6 kinds of therapy area, to make sure the treatment on the corner of eye, and other area on face.

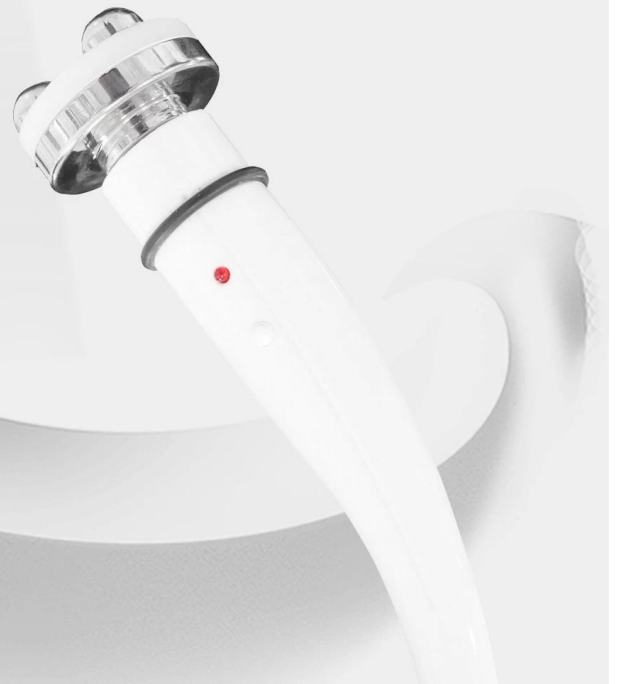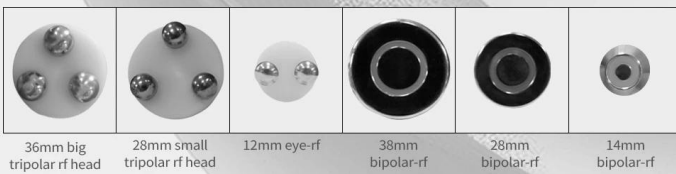

RF handle is suitable for 6 kinds of therapy area, to make sure the treatment on the corner of eye, and other area on face.

The RF handle is suitable for 6 kinds of therapy area, to make sure the treatment on the corner of eye, and other area on face.

The RF handle is suitable for 6 kinds of therapy area, to make sure the treatment on the corner of eye, and other area on face.

14 inch true color touch screen

Professional Screen, Easy operation



□□□□□□ □□□□ □□□□ □□□□ □□□□ 3-5 :□□□□□□ □□□□.

□□□□□□ □□□□: T/T (□□□□□□□□ □□□□□□□□ ,□□□□□□□□ □□□□□□)

□□□□□□ □□□□ □□□□ ,□□□□□□□□ □□□□□□□□ □□□□□□□□ :□□□□□□□□

□□□□□□□□ □□□□ □□□□ □□□□□□□□ :□□□□□□□□ □□□□, DHL/UPS/TNT

□□□□ □□□□ □□ □□□□□□□□□□□□ □□□□□□□□ :□□□□□□□□

Packing

safe and thoughtful packaging service



Built-in prevention Collision bubble

Two-year warranty for machine handle as consumable,
one-year warranty.

OEM&ODM service for distributor.

Tight packaging

Training: (Video+Manual+Online Service)
Can Keep You Operate It Easily.



3-5 business days delivery

Professional guide:
24 Hours Professional online customer service,