

More than 4-in-1

# 4IN1 MACHINE


Various combinations, as you wish





# Surprising Combination

Ultrasonic scalpel + gold microneedle + 980 +  
(choose one of 808, IPL, ND: YAG)

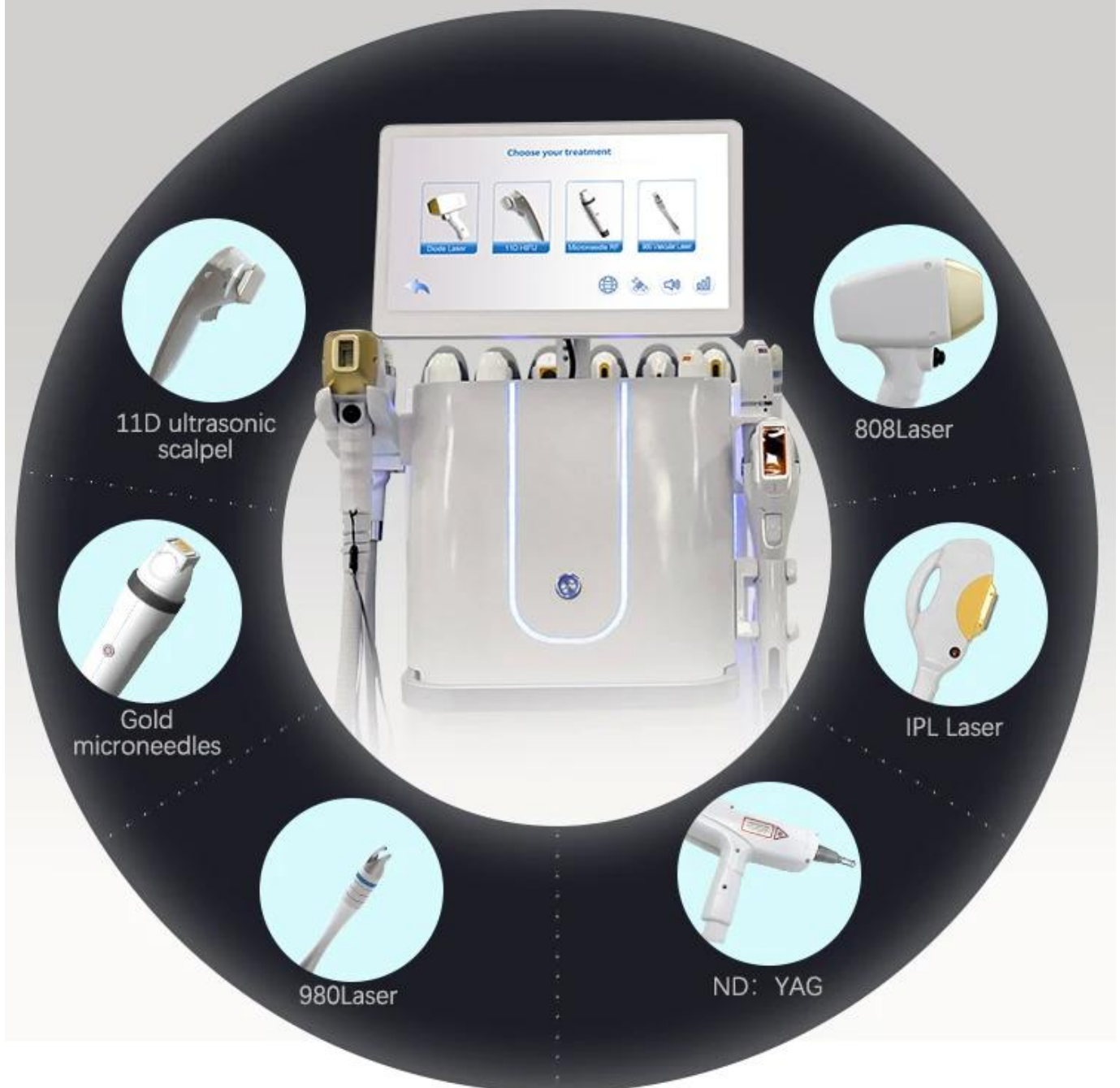

For this machine, totally has 4 handles. Match according to your needs.

1) Elight handle, use for:

Acne removal, vascular treatment, pigmentation therapy, remove freckles, skin rejuvenation, wrinkle removal, remove red blood silk, hair removal.

2) Diode Laser handle, use for:

Permanent hair removal.

3) ND YAG handle, use for:

Tattoo removal, tattoo on brow removal, lipline removal, skin rejuvenation, remove birth mark, otas nevus and dermal speckle.

4) HIFU handle, use for:

Wrinkle reduction contour shaping

5) Microneedle handle, use for:

Wrinkle reduction contour shaping□anti-wrinkle dissolve fat, lift and promote skin regeneration

6□980nm laser

Spider vein red bloodline removal, gray nail fungus treatment

## Product Description



# Various combinations

Meet your various needs



## Gold Microneedles

- Frequency 5 MHZ, better than most products on the market;
- Gold microneedles are better than non-insulated and ordinary ones;
- Four types of treatment heads:

## 980 Laser

- Using a new generation of imported lasers, the treatment effect is significantly improved.



## IPL Laser

- High-quality IPL: handpiece with super cooling system
- Optional Spot Size: 15\*50mm, 15\*35mm, 12\*30mm
- 7 filters: 420/510/560/585/640/690/755nm  
Standard: 420/510/560/690nm

## ND YAG

- 5 treatment heads:  
A solution to a variety of skin problems



## About Handle

# DIODE LASER

Suit to any color of skin

- **PAINLESS:** 755/808/1064nm diode laser for all type skin hair removal
- **FASTER:** Large area spot sliding operation
- **NO TRAUMA:** There is no trauma to skin tissue
- **SAFETY:** To avoid the light of the deputy to reduce the difficulty of operation



755nm



808nm



1064nm

Diode laser suits *for all skin* and hair types treatment.

Imported Laser Bar. 50 million shots time.



# FAC+Sapphire lenses:

Energy transfer efficiency up to 99.9%



**Without FAC  
The Beam Diverges Rapidly**

Ata working distance of 32mm  
The energy density about  
2J/cm<sup>2</sup> @ 20ms



**Add FAC: Beam Focus**

Areslite Ata working distance of 32mm .  
The energy density is 10J/cm<sup>2</sup>@ 20ms



# CLINICAL EFFECT



# Elight Handle Design



## Specification

**Optional Spot Size:**  
15\*50mm

**Imported UK lamp**

**Long shots time:**  
150,000--300,000 times

**7filters:**  
420/510/560/585/640/690/755nm

**Super cooling system**



Elight for acne removal, vascular treatment, pigmentation therapy, remove freckles, skin rejuvenation, wrinkle removal, remove red blood silk, hair removal.

|   |   |   |   |  |   |   |
|---|---|---|---|--|---|---|
|  |  |  |  |  |  |  |
| <b>420nm</b><br>Acne removal<br>Vascular treatment                                | <b>510nm</b><br>Pigmentation therapy, Remove freckles                             | <b>560nm</b><br>Skin rejuvenation<br>wrinkle removal                              | <b>585nm</b><br>Remove Red blood silk   | <b>640nm</b><br>Hair removal (light skin)  | <b>690nm</b><br>Hair removal (brown skin)   | <b>755nm</b><br>Hair removal (black skin)   |
|  |  |  |  |  |  |  |
|  |  |  |  |  |  |  |

# 5 therapeutic **Heads**

■ 532 & 1064: removing tattoo, eyeliner and lip line

■ 1320: black-faced doll whitening and whitening skin

■ Adjustable 532 & 1064: removing large-area tattoos



Adjustable  
532nm



Adjustable  
1064nm



Black doll  
1320nm



Fixed  
532nm



Fixed  
1064nm

ND YAG handle, use for tattoo removal, tattoo on brow removal, lipline removal, skin rejuvenation, remove birth mark, otas nevus and dermal speckle.

# Working Principle





# 11 Heads!

## 7 Cutter Heads, 4 Cannon Heads



### SIZE OF CARTRIDGE

1.5mm,2.0mm,3.0mm,4.5mm for face treatment  
6.0mm.9mm.13mm for body treatment.

11D Cartridges L7-1.5mm,MF2.0mm,L7-3.0mm,L4-4.5mm,MF6.0mm,  
MF9.0mm,MF13.0mm,DL1.5,DL2.0,DL3.0,DL4.5


1. Remove wrinkles on around forehead, eyes, mouth, etc.
2. Lifting and tightening both cheeks skin.
3. Improving skin elasticity and shaping contour.
4. Improving jaw line, reducing "marionette lines".
5. Tightening the skin tissue on forehead, lifting the eyebrows lines.
6. Improving skin complexion, making the skin delicate and bright.
7. Match with injecting beauty like Hyaluronic acid, collagen, to resolve more aging problem.
8. Removing neck wrinkles, protecting neck aging.

## ***Comparison before and after use***






# Handle details

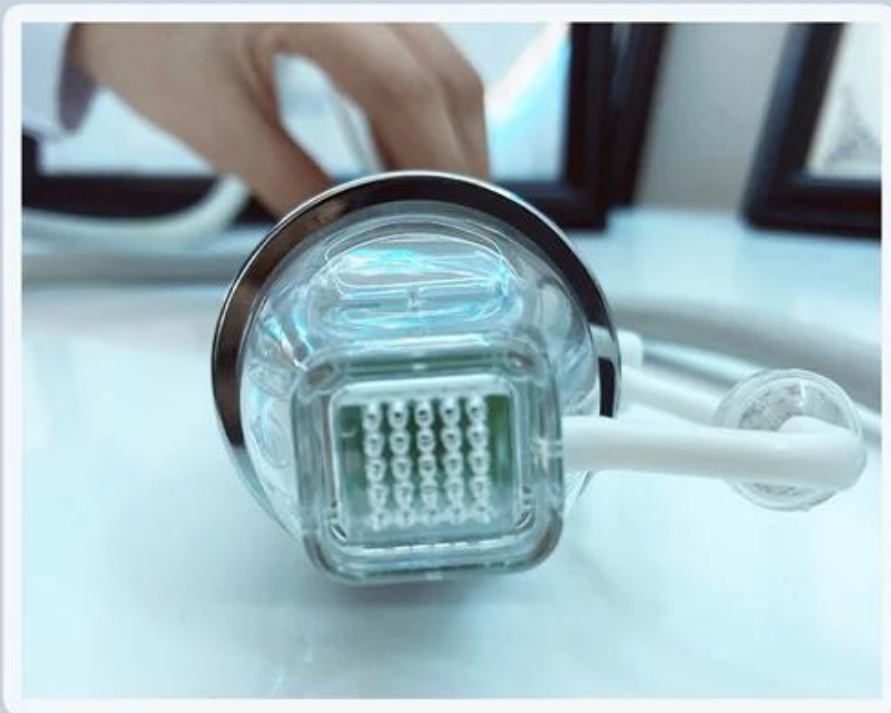


1. Non insulated needles: equal therapy to epidermis and dermis layer 2. Stepping motor type: needle inserts to skin smoothly without shock 3. Gold plated needles: high biocompatibility, suits metal allergy patient \*. Anti-wrinkle, firming the skin, improve pseudo wrinkles, dissolving fat, shaping and lifting. \*. Rapidly brighten the skin and make it more tender. \*. Actively promote facial lymph circulation, solve the problem of skin edema. \*. Lifting and tightening the skin, effectively solving repairing stretch marks. \*. Removing black rim of the eye, eye bags and wrinkles around the eye. Shrinking the pores, repairing blain scar, calming skin.

# WORKING PRINCIPLE

Make Your Beauty Business Easier

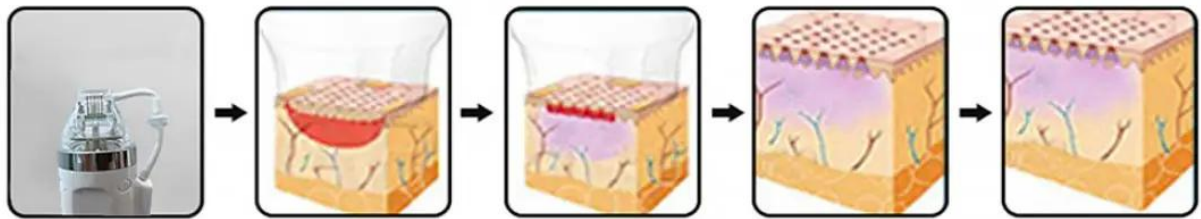

Microneedle delivers thermal energy to the treatment area more safely and effectively allowing treatment regardless of skin type. The laser energy is mostly absorbed by melanin, especially in darker skin which contains a higher melanin content.

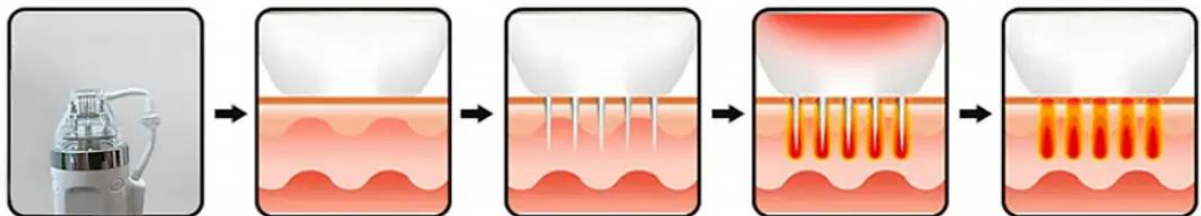

For Micro Needle Microneedle treatments, Effective, direct thermal stimulation to remodel collagen with less epidermis damage. Micro Microneedle penetrates target area from epidermis to dermis without energy loss; stimulates and efficiently remodels disordered collagen.

# PROFESSIONAL AND FAST

CHANGE ➤



980nm laser beam is focused onto a 0.5-3mm diameter range, in order to enable more focused energy to reach the target tissue, while avoiding burn the surrounding skin tissue. The transcutaneous irradiation (with a penetration of 1 to 2 mm in the tissue) causes tissue selective absorption by hemoglobin (hemoglobin is the main target of the laser) Spider vein red bloodline removal, gray nail fungus treatment

# TREATMENT EFFECT T

BEFORE & AFTER



OPERATION SHOWS





# 14 inch 4K true color touch screen

Operation interface, independently designed to ensure  
safety and easy operation



1. 14inch big touch screen, humanized designed

2. Male/femal, body parts, 6 skin types for option.

3. English+your native language.

4. Free OEM program service(logo/system...)

If you have any questions about beauty equipment, looking for partners and beauty entrepreneurship, please consult MLKJ!

**Tel/Whatsapp: +86 15965364558**

**E-mail: [sale@mingliangkj.com](mailto:sale@mingliangkj.com)**