

980+ microneedle multi-function machine

- ✓ Vascular lesions
- ✓ Proliferative lesion
- ✓ Age spots and sun spots
- ✓ Spider veins
- ✓ Face lift
- ✓ Skin rejuvenation
- ✓ Reduce wrinkles
- ✓ Scar remover
- ✓ Skin tightening
- ✓ Acne treatment




980nm Diode Laser + Micro Needle Applications Therapeutic range

- 1) Vascular lesions □ Proliferative lesion
- 2) Spider vein removal
Telangiectasia and blood vessels treatment
- 3) Face lift □ skin rejuvenation
- 4) Acne treatment □ reduce wrinkles □ scar remover □ skin tightening
- 5) Remove speckle, age spots, sunburn, pigmentation

EXFU

Focus on repairing skin problems

Clear effect, strong pertinence.

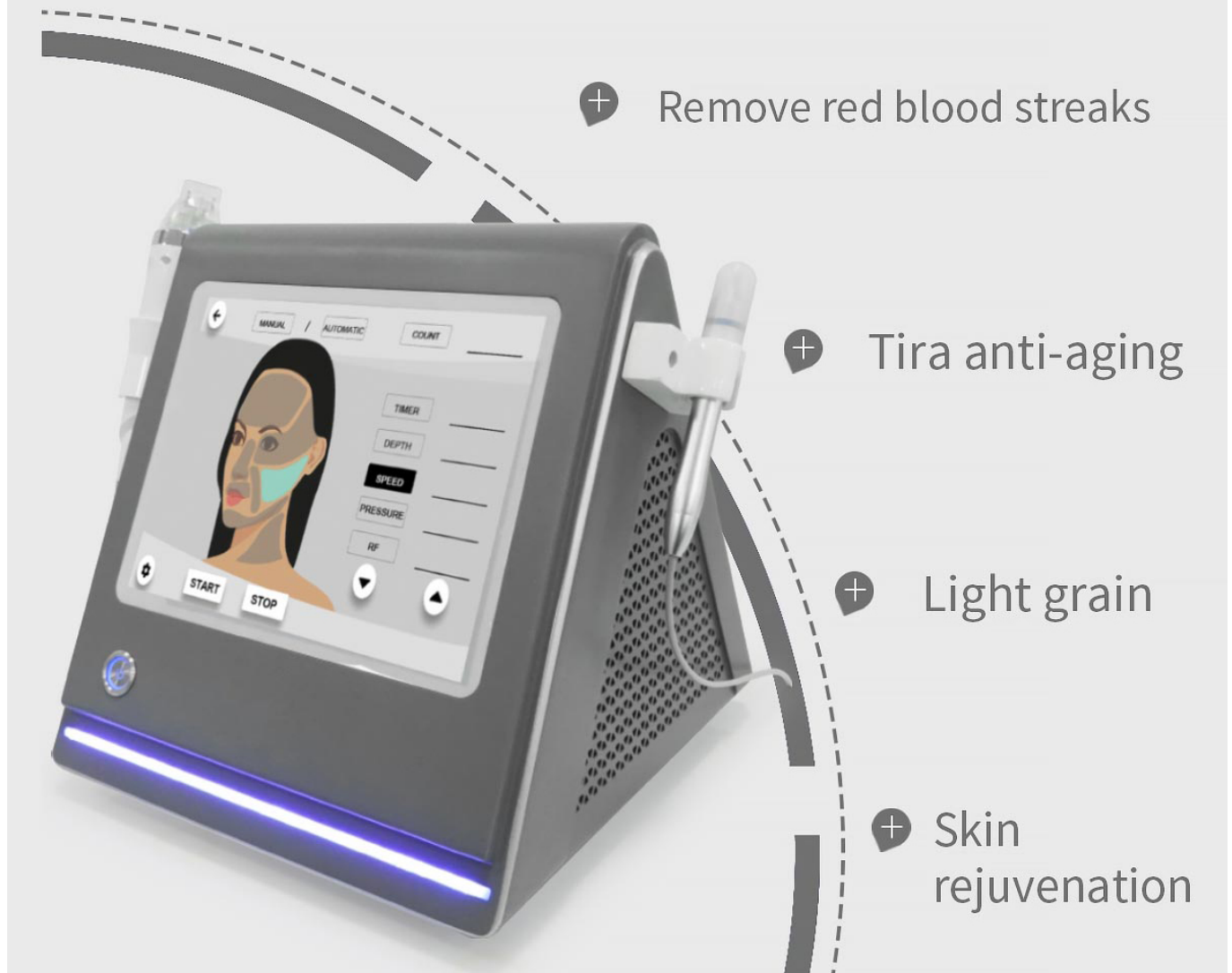

2 in 1 980nm Diode Laser + Micro Needle

1. 980nm laser beam is focused onto a 0.2-0.5mm diameter range, in order to enable more focused energy to reach the target tissue, while avoiding burn the surrounding skin tissue. The transcutaneous irradiation (with a penetration of 1 to 2 mm in the tissue) causes tissue selective absorption by hemoglobin (hemoglobin is the main target of the laser).

2. RF Microneedling treatment uses ultrafine needles and radiofrequency (RF) energy to penetrate the top layer of the skin. This allows for the regeneration of elastin and collagen to the treated areas, resulting in a revitalized appearance of skin.

Microneedle 980 laser





RF Microneedling Handle Display


Handle details




The Exfu RF Microneedling machine can tighten skin, and ensure firmness, through soft-tissue coagulation. This allows for the regeneration of elastin and collagen to the treated areas, resulting in a revitalized appearance of skin.

WORKING PRINCIPLE

Make Your Beauty Business Easier

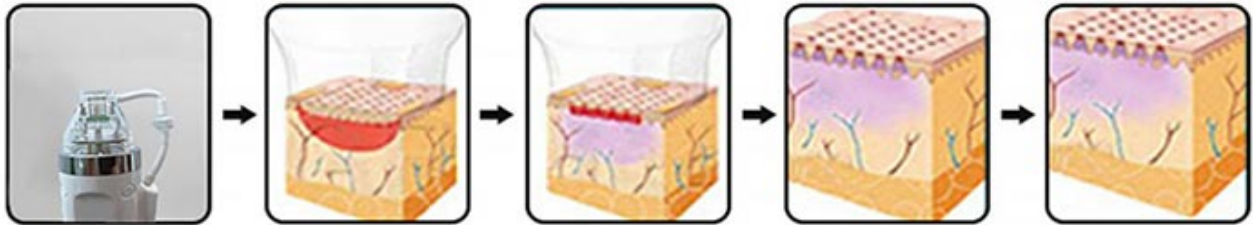


Microneedle delivers thermal energy to the treatment area more safely and effectively allowing treatment regardless of skin type. The laser energy is mostly absorbed by melanin, especially in darker skin which contains a higher melanin content.

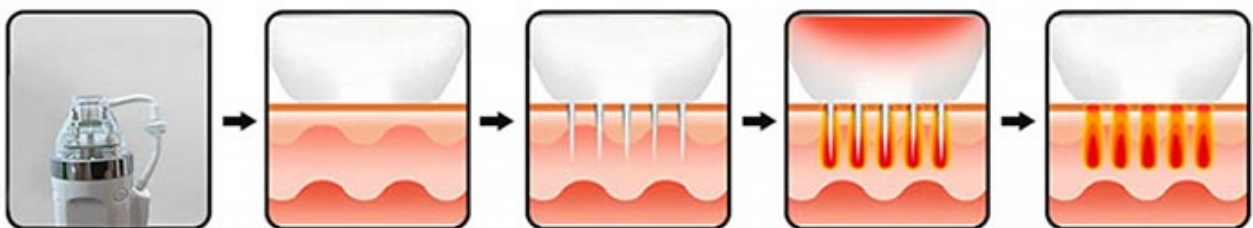

For Micro Needle Microneedle treatments, Effective, direct thermal stimulation to remodel collagen with less epidermis damage. Micro Microneedle penetrates target area from epidermis to dermis without energy loss; stimulates and efficiently remodels disordered collagen.

How does treat for skin revitalization? Wake up your skin's natural collagen & elastin production. Frequency microneedling treatment uses ultrafine needles and radiofrequency (RF) energy to penetrate the top layer of the skin and trigger the body's natural healing process to regenerate new collagen and elastin. Unlike some other skin revitalization treatments, these treatments can be performed on all skin types, anywhere on the body and anytime of the year.

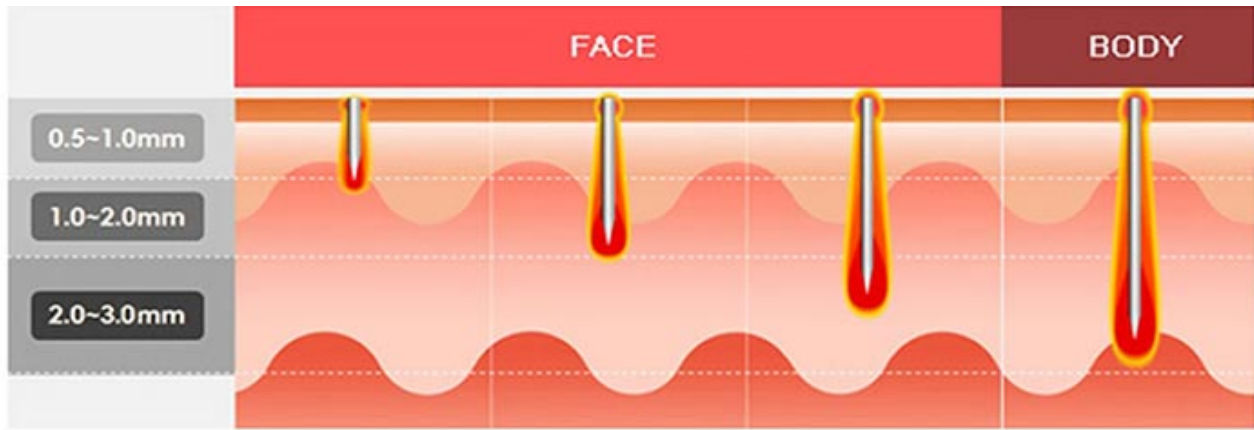

Control the treatment depth from 0.5-3mm by screen

TREATMENT EFFECT

Make Your Beauty Business Easier



What are the treatment areas?

Both face and body are able to be treated. The treatable areas are as follows: Full Face/Forehead/Crow's feet/Upper Cheek/Lower Eyes/Lower Cheek/Nose/Upper & Lower Lip/Chin/Jawline/Neck/Decolletage/Elbows/Abdomen

980nm Diode Laser Handle display

PROFESSIONAL AND FAST

CHANGE 



- **Laser power 15w**

Our 980 machine adopts 15w low power, adjustable power and 30w application in medical devices.



- **Adjustable light spot**

Adjusting ring with four adjusting light spots, the thicker the adjusting ring, the larger the light spot, and the thinner the adjusting ring, the smaller the light spot.



- **Life 10,000 hours**

The service life is up to 10,000 hours, which can be used for 24 hours continuously.



- **No trauma**

Max directional conduction patent technology is adopted to act on the capillary in a targeted manner, breaking up hemoglobin in the capillary in an instant, shrinking and solidifying the expanded capillary in an instant, and the red blood filaments disappear.

- **Two therapeutic heads**


The two treatment heads are respectively used for treating red blood silk and onychomycosis.




Spider Vein Removal 980 Diode Vascular Laser Machine
Diode Laser 980nm repair and remove gray nails

WORKING PRINCIPLE

The pulse width can be used to warm up the blood vessels so that they can be heated and condense to avoid the purple

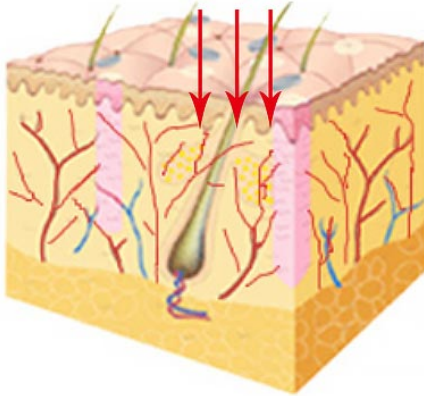

1 The hemoglobin in the capillary rupture of the moment

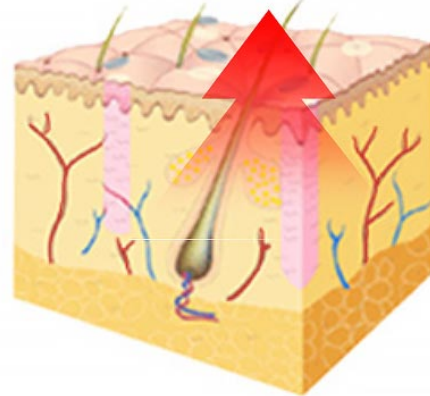

2 The hemoglobin was removed after rupture

TREATMENT EFFECT

The pulse width can be used to warm up the blood vessels so that they can be heated and condense to avoid the purple

OPERATION SHOWS




BEFORE & AFTER



TREATMENT EFFECT T

BEFORE & AFTER

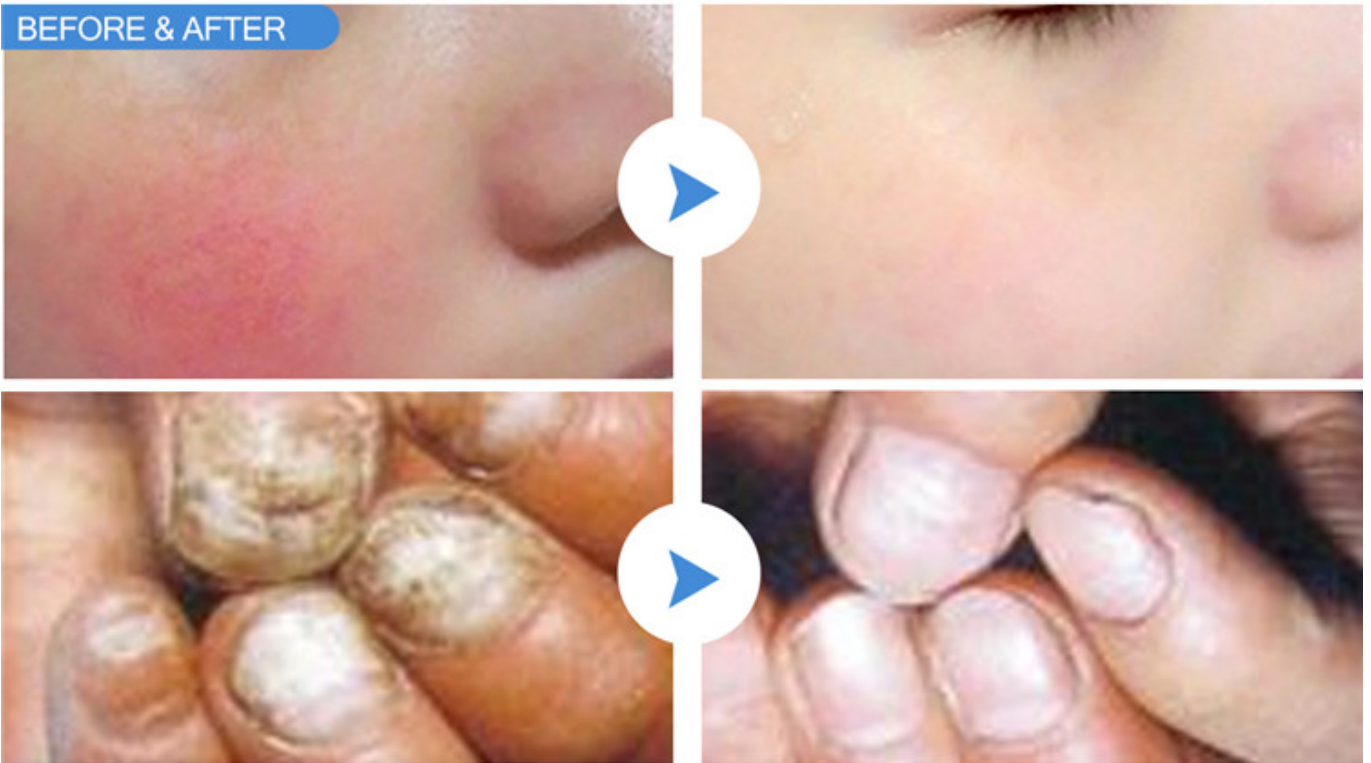

OPERATION SHOWS



Machine features:

- * No consumable parts, the machine can work 24 hours a day.
- * Not damage the epidermis at all. * Painless, no-invasive, no downtime.
- * Portable design, easy for transport.
- * Top equipped spare, 980nm laser with 4 rings to adjust spot size+gold micro needle machine with 4 heads for skin resurfacing.
- * Multiple languages can be selected, we can make the free logo according to your requirements.



One-key power-on design

Screen Display

4K SCREEN DISPLAY

Make Your Beauty Business Easier



1. 14 inch big touch screen, humanized designed
2. Male/femal, body parts, 6 skin types for option.
3. English+your native language.
4. Free OEM program service(logo/system...)

Certificate display

Certificate is the best guarantee for you





Company Profile



15 Years
Laser Beauty Device
Manufacturer

MLKJ is a 15-year company integrating beauty instrument design, development, production and sales. Real factory, real production process, real evaluation! Welcome to our company or video visit! Only when you choose a real partner can the product quality be guaranteed!

Our Showroom



Production Process



Warehouse spare parts



Parts assembly



Precision machining



Handle fabrication



Product assembly process



Quality inspection
of finished products



Product packaging



Packing

Safe and thoughtful packaging service



Tight packaging

The machine shall be made of waterproof materials and packaged

**Anti collision foam board
inside aluminum box**



3-5 business days delivery

Professional guide:

24 Hours Professional online customer service