


Short pulse laser: duration \leq nanoseconds (10e-6 seconds) Best laser for tattoo removal

Features:

- \square Using the mechanical sound pressure blasting effect of light, a single emission duration is short, high energy, peak power;
- ☐ Pulse width is generally not adjustable;
- ☐ Short pulse is realized by Q-switching technology

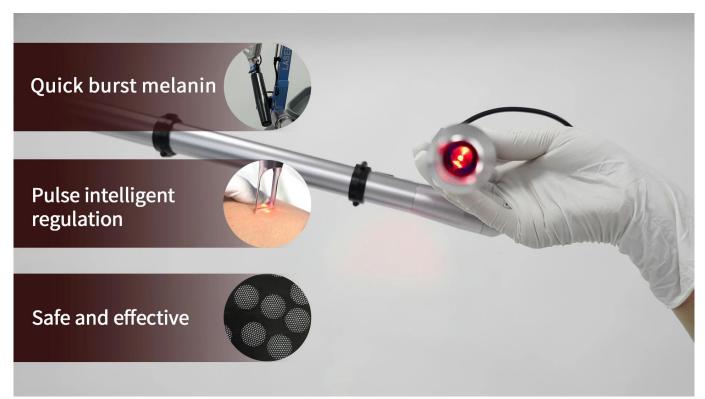


Product Information




With 7 guiding light arms, the light is more symmetrical! Energy density (J/cm2) is the decisive parameter of Q-switched laser therapy. Q-switched laser usually can not adjust the pulse width, set the same light energy value (mJ), if the smaller the area of action (the smaller the spot diameter), the greater the energy density, the more intense the skin reaction, which is the meaning of "focusing"

About Handle

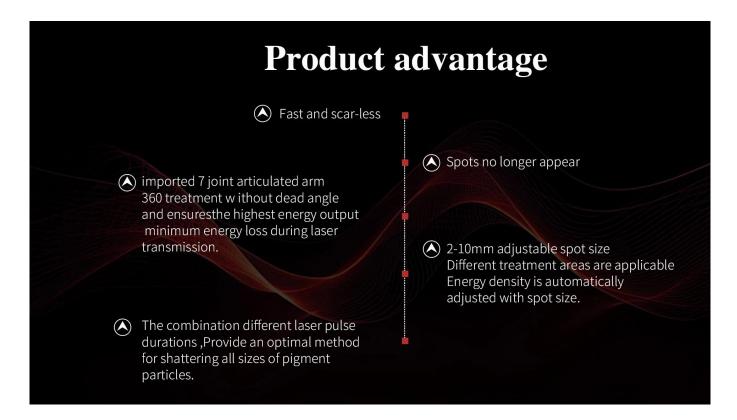

Treatment Effect

Features: 1. High performance: The Ruby laser system has a high energy output and high power density, enabling quick and effective surgery and treatment. 2. Safety: The system uses advanced laser technology to accurately control the depth and intensity of the laser beam, minimizing damage to the surrounding tissue. 3. Multi-functional: Ruby laser system can be used for a variety of surgeries and treatments, including skin removal, vascular disease treatment, scar repair, etc. 4. High precision: The system has high precision laser beam positioning and control functions, can accurately handle the target tissue, improve the success rate of surgery and treatment effect. 5. Portability: Ruby Laser system adopts compact design, easy to carry and move, suitable for different clinical environments.

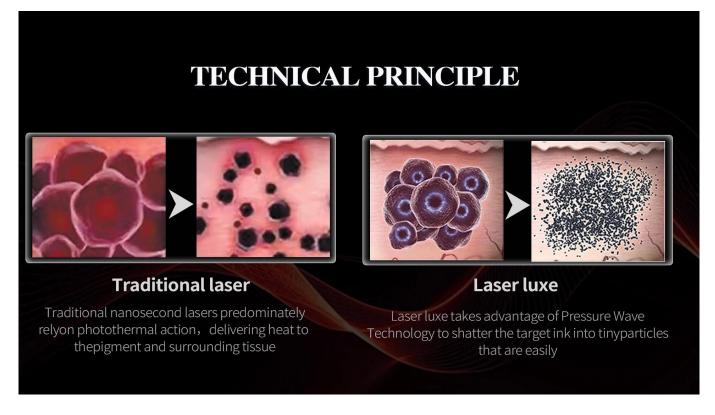

755 waveleng this used for whitening and yellowing quickly remove dust pigments, and easily restore your white skin The 532 nm picosecond laser has the best effect on yellow tattoos, followed by red tattoos. The 1064 nm picosecond laser works best for black tattoos. 532 Wavelength is used to treat epidermal spots such as freckles, coffee spots, photo-aging spots, acne marks and wrinkles, as well as red and yellow tattoos. 1064 used to treat dermal spots such as brown-

cyanmoles,Otamoles,chloasma,pigmenta-tion,whitening and brightening to remove yellow,black and blue tattoos can be washed,used to treat acne,reduce pores,and tightenskin,To achieve the effect of skin rejuvenation.




Parameters Tal	ble	
Device Name	Laser luxe	
Wavelength	1064nm 755nm 532nm	
Screen Size	12 inch 4K Style touch screen	
Pulse Duration	\sim 7ns@1064nm \sim 6.5ns@532nm	
Frequency	1-10hz	
Spot size	2-10mm adjustable	
Energy	1064/1000mj 755/1000mj 532/1000mj	
Beam of light transmission	7 articular-arm of light guide transmitting	
Indicator of aiming light	650 red diode laser	
Cooling	Water cooling+Air cooling	
Voltage	AC220V/110V 50HZ/60HZ	
Demension of machine	93.6*92.7*39cm	
Package	Aluminum box	
Warranty	2 years	

Skin Whitening Picosecond Laser Q-swithced Nd Yag Laser Tattoo Removal Machine's Accessories: In addition, we also have various accessories in order to give our customers excellent service. Laser Q-swithced nanosecond short pulse 755 532 1064 vascular therapy scar repair picosecond laser tattoo removal machine





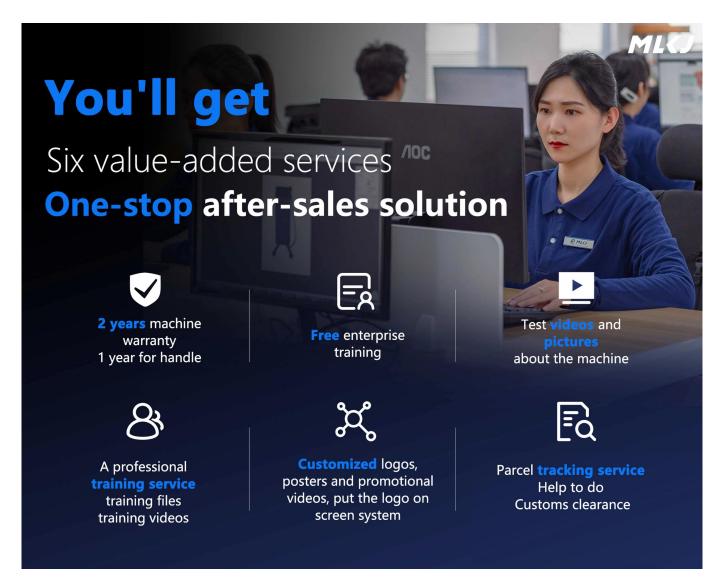

After Service:

- *Considering the long distance between us with our clients, and in order to serve our clients much better, *Our company build a professional technology after service team.
- *Our technician will provide timely solution for any machine problem and prepare the documents for you within 24 hours.
- *You can also contact us by phone, online etc..

Learn more product information and beauty knowledge, welcome to consult us! 7*24 hours online service!